
# TERUMO

# SURFLO™

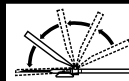**WINGED INFUSION SET WITH  
NEEDLE PROTECTION**

**CE**  
**MD 0197**

en	WINGED INFUSION SET WITH NEEDLE PROTECTION	2-3
fr	MICROPERFUSEUR À AILETTES AVEC PROTECTEUR D'AIGUILLE	4-5
de	PERFUSIONSBESTECK MIT FLÜGELKANÜLE UND NADELSCHUTZ	6-7
es	PALOMILLA DE INFUSIÓN CON PROTECCIÓN DE LA AGUJA	8-9
pt	DISPOSITIVO COM ASAS PARA INFUSÃO COM PROTETOR DE AGULHA	10-11
it	AGO A FARFALLA CON COPRIAGO DI PROTEZIONE	12-13
nl	INFUSIESET MET VLEUGELNAALD EN NAALDPROTECTOR	14-15
sv	INFUSIONSSET MED VINGAR OCH KANYLSKYDD	16-17
da	INFUSIONSSÆT MED VINGE OG KANYLEBESKYTTELSE	18-19
no	VINGET INFUSJONSSETT MED NÅLEBESKYTTELSE	20-21
fi	SIIPINEULA INFUUSIOSETTI NEULASUOJUKESELLA	22-23
el	ΣΕΤ ΕΓΧΥΣΗΣ ΜΕ ΠΕΤΑΛΟΥΔΑ ΜΕ ΠΡΟΣΤΑΤΕΥΤΙΚΟ ΒΕΛΟΝΑΣ	24-25
ru	УСТРОЙСТВО ДЛЯ ИНФУЗИИ С ИГЛОЙ «БАБОЧКА» И ПРОТЕКТОРОМ ИГЛЫ	26-27
pl	MOTYLKOWY ZESTAW INFUZYJNY Z ZABEZPIECZENIEM IGLY	28-29
hu	SZÁRNYAS INFÚZIÓS KÉSZLET TŰVÉDELEMMEL	30-31
cs	KRIDÉLKOVÝ INFUZNÍ SET S OCHRANOU JEHLY	32-33
sk	KRÍDIELKOVÝ INFÚZNY SET S CHRÁNIČOM IHLY	34-35
tr	İĞNE KORUMALI KANATLI INFÜZYON SETİ	36-37
et	LIBLIKÕELA INFUSIOONIKOMPLEKT KOOS NÕELAKAITSEGA	38-39
lv	TAURIŅVEIDA INFŪZIJAS KOMPLEKTS AR ADATAS AIZSARDZĪBU	40-41
lt	PETELIŠKĖS FORMOS INFUZIJOS RINKINYS SU ADATOS APSAUGA	42-43
sl	INFUZIJSKI SET S KRILCI TER ZAŠČITO IGLE	44-45
sr	INFUZIONI SET SA KRILCIMA I ZAŠTITOM ZA IGLU	46-47
ro	MICROPERFUZOR CU ARIPIOARE ŞI PROTECTOR DE AC	48-49
bg	ИНФУЗИОНЕН НАБОР С КРИЛЦА И ПРЕДПАЗИТЕЛ НА ИГЛАТА	50-51
uk	ІНФУЗІЙНА СИСТЕМА-МЕТЕЛИК ІЗ ЗАХИСТОМ ГОЛКИ	52-53
hr	KOMPLET ZA INFUZIJU S KRILCIMA I ZAŠTITOM IGLE	54-55


# Instructions for Use

**ENGLISH**
**SYMBOL EXPLANATION**

<b>LOT</b>	Batch code	Do not re-use	<b>MD</b>	Medical Device
<b>REF</b>	Catalogue number	Do not use if package is damaged and consult instructions for use	Non pyrogenic	
<b>Rx only</b>	Caution: Federal law (USA) restricts this device to sale by or on the order of a physician	Fragile, handle with care	Single Sterile Barrier System	
	Consult Instructions for Use	Keep away from sunlight	Sterilized using ethylene oxide	
	Contents	Keep dry	Temperature limit	
	Date of Manufacture	Manufacturer	Use-by date	
	Do not resterilize			

**PRODUCT DESCRIPTION**

Non toxic

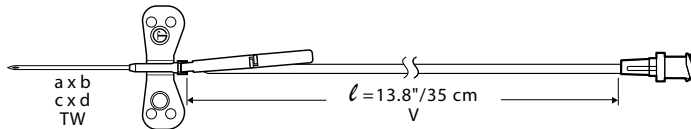Catalogue number	Needle Gauge (a)	Needle Length (b)	Needle Outer Diameter (c)	Needle Length (d)	Wall Thickness	Tube Length (l)		Dead Space Volume (V)
SV-S23WL35	23G	3/4"	0.6 mm	19 mm	Thin Wall (TW)	13.8"	35 cm	≤ 0.20 ml
SV-S23BL35			≤ 0.40 ml					
SV-S25WL35	25G	3/4"	0.5 mm	19 mm	Thin Wall (TW)	13.8"	35 cm	≤ 0.20 ml
SV-S25BL35			≤ 0.40 ml					

**INTENDED PURPOSE**

The Surflo™ Winged Infusion Set with Needle Protection (Surshield™) is intended to access the peripheral vascular system, for single-dose or short-term intravenous administration of fluids and/or withdrawal of blood specimen using a syringe, a luer adapter or another compatible/appropriate device. Additionally, after withdrawal of the needle from the patient's vein, the safety shield (needle protection (Surshield™)) shall be manually activated to cover the needle to minimize risk of accidental needle stick.

**INDICATIONS**

The Surflo™ Winged Infusion Set with Needle Protection (Surshield™) is for general application - for treatment (administration of fluids) or for diagnosis (withdrawal of blood specimen).

**CONTRA-INDICATIONS**

No contra-indications.

**PATIENT TARGET GROUP**

Intended for general application.

**INTENDED USERS**

Healthcare professional or lay person.

**CLINICAL BENEFIT**

The Surflo™ Winged Infusion Set with Needle Protection (Surshield™) has an indirect clinical benefit (indirect medical effect) since it is used for administration of fluids or for withdrawal of blood specimen. For general application. The safety shield (needle protection (Surshield™)) has an indirect clinical benefit since it prevents needle stick injury.

**WARNINGS**

- Do not use if unit package is damaged.
- Use immediately after opening the unit package.
- The needle is made of stainless steel containing nickel and cobalt. Cobalt is classified as CMR\* 1B and is present in a concentration above 0.1% weight by weight. Current scientific evidence supports that medical devices manufactured from stainless steel alloys containing cobalt do not cause an increased risk of cancer or adverse reproductive effects.
- \* CMR = Carcinogenic, mutagenic or toxic to reproduction (CLP Regulation EU 1272/2008)

**PRECAUTIONS**

- After use, secure the needle into the safety shield according to the directions for use.
- If the needle is bent or damaged, no attempt should be made to straighten the needle or use the product.
- After use, dispose of safely as medical waste in a sharps disposal container and/or according to health institution policies. The product is biohazardous and is physically hazardous due to its sharp edge.

**DIRECTIONS FOR USE**

Aseptic technique, proper skin preparation and continued protection of the site are essential. Observe Universal Precautions on ALL patients.

**PRECAUTION** Keep hands behind the needle at all times during use and disposal.

**SET UP FOR ADMINISTRATION**

1. Open the package.
2. a. Connect a syringe or another compatible/appropriate device and secure.
  - b. Proceed with priming and confirm that administration fluid comes from the needle.

**SET UP FOR BLOOD COLLECTION**

1. Open the package.
2. Connect a syringe, a luer adapter or another compatible/appropriate device and secure.

**VENIPUNCTURE FOR ADMINISTRATION / BLOOD COLLECTION**

3. Flip the safety shield back, away from the needle towards the tubing. Grasp wings.
4. Remove the needle protector.
 

**Caution:** Care should be taken not to touch the needle.
5. Perform venipuncture following facility protocol and confirm proper positioning in the vessel (*Fig.1*).
6. Secure the position of winged infusion set per facility protocol.

**AFTER USE**

7. Remove tape, if present, from wings.
8. Flip the safety shield forward towards the needle (*Fig. 2*). Grasp the wing with your thumb on top of the shield and your index finger under the wing. Completely remove the needle from the puncture site and apply digital pressure to the site using a sterile gauze pad (*Fig. 3*).
9. While maintaining the same grip on the device (per step 8 above), activate the safety measure by performing either of the following methods:

**► Recommended (one-handed) method:**

(Reference OSHA 1910.1030 (d)(2)(vii)(B))







Slide your thumb forward over the shield while pinching your finger and thumb together (or press the safety shield against a hard surface such as a bedside table) until an audible click is heard and the safety feature is visually confirmed to be activated (*Fig. 4a*).

**OR**
**► Alternate (two-handed) method:**

While keeping both hands behind the needle, grasp one of the wing tabs with your free hand to stabilize the device. (If necessary, reposition your original grip per step 8 being careful not to press against an exposed area of the needle.) Then, slide your thumb forward over the shield while pinching your finger and thumb together until an audible click is heard and the safety feature is visually confirmed to be activated (*Fig. 4b*).

**Caution:** Needle must be carefully and firmly placed in the lock position of the safety shield.

10. Do not attempt to deactivate the safety device by separating the needle from the safety shield.
11. Dispose of safely as medical waste according to health institution policies (*Fig. 5*).

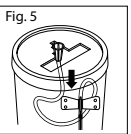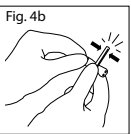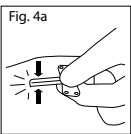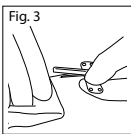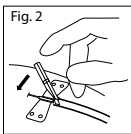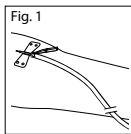**PRECAUTIONS FOR STORAGE**

Store between 2°C and 30°C. Keep dry. Keep away from sunlight. Fragile, handle with care.

**REPORT OF INCIDENT**

If during the use of this device or as a result of its use, a serious incident has occurred, please report it to the manufacturer and to your national authority.

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