

Senri®

0.018" RX PTA Balloon Dilatation Catheter

Low Profile PTA Balloon combining Finesse and Pushability



Low profile semi-compliant balloon, hydrophilic coated

- ▶ Featuring an extremely low crossing profile, by a 3-folded ($\varnothing 3.0-6.0\text{mm}$) and 4-folded ($\varnothing 7.0-8.0\text{mm}$) balloon with high flexibility for successful crossability and superior trackability.
- ▶ The distal 45 cm catheter segment is hydrophilic coated for swift and easy navigation, even through complex anatomy.

Dedicated shaft design

- ▶ Combines the pushability of a stainless steel hypotube shaft with the kink resistance of a polymer shaft, even for cross-over procedures in sharp angulated bifurcations.
- ▶ Maintains very fast deflation times.

Broad range of sizes

- ▶ Balloon diameters from 2-8mm and lengths up to 150mm - suitable for a great variety of peripheral interventional procedures.


Senri®

0.018" RX PTA Balloon Dilatation Catheter

Specifications

Balloon material: Polyamide
 Balloon folding: ▶ 3-folded (ø 3.0 - 6.0 mm)
 ▶ 4-folded (ø 7.0 - 8.0 mm)

Maximum guide wire: 0.018" (0.46 mm)
 Effective catheter length: 150 cm

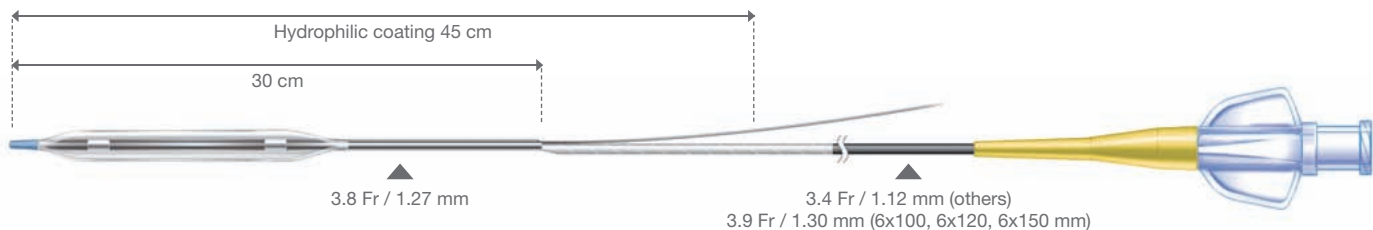

Order details

Balloon diameter	Code number						Minimum sheath size	Minimum guiding catheter	Nominal pressure	Rated burst pressure
	Balloon length									
	40 mm	60 mm	80 mm	100 mm	120 mm	150 mm				
2.0 mm	BD-S2040L	BD-S2060L	BD-S2080L	BD-S20100L	BD-S20120L	BD-S20150L	4 Fr (1.33 mm)	6 Fr (2.0 mm)	8 atm/ 811 kPa	14 atm/ 1 419 kPa
3.0 mm	BD-S3040L	BD-S3060L	BD-S3080L	BD-S30100L	BD-S30120L	BD-S30150L				
4.0 mm	BD-S4040L	BD-S4060L	BD-S4080L	BD-S40100L	BD-S40120L	BD-S40150L				
5.0 mm	BD-S5040L	BD-S5060L	BD-S5080L	BD-S50100L	BD-S50120L	BD-S50150L				
6.0 mm	BD-S6040L	-	-	-	-	-	5 Fr (1.67 mm)	7 Fr (2.33 mm)		
	-	BD-S6060L	BD-S6080L	BD-S60100LH	BD-S60120LH	BD-S60150LH				
7.0 mm	BD-S7040L	BD-S7060L	BD-S7080L	-	-	-	6 Fr (2.0 mm)	8 Fr (2.67 mm)		
8.0 mm	BD-S8040L	BD-S8060L	-	-	-	-				



Terumo Corporation
 44-1, 2-chome, Hatagaya
 Shibuya-ku, Tokyo 151-0072,
 Japan
 Tel.: +81 3 3374 8111
 Fax: +81 3 3374 8196

European Headquarters

Terumo Europe N.V.
 Interleuvenlaan 40
 3001 Leuven, Belgium
 Tel.: +32 16 38 12 11
 Fax: +32 16 40 02 49

**Terumo Europe
 Interventional Systems
 Liaison Office**
 Immeuble Krebs
 8, rue Jean-Pierre Timbaud
 CS 10720
 78067 Saint Quentin
 en Yvelines Cedex,
 France
 Tel.: +33 139 30 95 75
 Fax: +33 139 30 95 70

European Offices

**Terumo Europe N.V.
 Benelux Sales Division**
 Interleuvenlaan 40
 3001 Leuven,
 Belgium
 Tel.: +32 16 39 25 80
 Fax: +32 16 39 25 99
 The Netherlands:
 Tel.: 0800 0220396
 Fax: 0800 0220414

Laboratoires Terumo France S.A.
 Parc Ariane - Bât. Uranus
 1, rue Hélène Boucher
 78284 Guyancourt Cedex,
 France
 Tel.: +33 130 96 13 00
 Fax: +33 130 43 60 85

Terumo Italia S.r.l.
 Via Simone Martini, 143 / 145
 00142 Roma RM, Italy
 Tel.: +39 0651 96 14 20
 Fax: +39 0650 30 407

Terumo UK Ltd
 Tamesis, The Causeway
 Egham, Surrey
 TW20 9AW, United Kingdom
 Tel.: +44 1784 476 200
 Fax: +44 1784 476 201

Terumo Deutschland GmbH
 Ludwig-Erhard-Strasse 6
 65760 Eschborn, Germany
 Tel.: +49 6196 80 230
 Fax: +49 6196 80 23 200

Terumo Deutschland GmbH
 Zweigniederlassung Spreitenbach
 Bodenackerstrasse 3
 8957 Spreitenbach, Switzerland
 Tel.: +41 56 419 10 10
 Fax: +41 56 419 10 11

Terumo Europe España SL
 Avda. Juan Carlos I, Nº 13 - 7ª Planta
 Edificio Torre La Garena
 28806 Alcalá de Henares (Madrid), Spain
 Tel.: +34 9021 01 298
 Fax: +34 9021 01 358

Terumo Sweden AB
 Nya Varvet, Byggnad 90
 Sven Källfelts Gata 18
 426 71 Västra Frölunda, Sweden
 Tel.: +46 3174 85 880
 Fax: +46 3174 85 890

**Terumo Denmark
 Filial of Terumo Sweden AB**
 Tel.: +45 7020 93 80
 Fax: +45 7020 94 80

**Terumo Europe N.V.
 Moscow Representative Office**
 BC Northern Tower
 13Fr., 10 Testovskaya Street
 Moscow 123317, Russia
 Tel.: +7 495 988 4740
 Fax: +7 495 988 4739

**Addresses subject to change,
 please consult our Website.**

© Registered Trademark
 Published by Terumo Europe N.V.