



At Terumo Interventional Systems, we constantly **work to refine and perfect our products** so that interventionalists can do more. That is why we support great thinking that pushes back the boundaries of our field.

# PUSHING BOUNDARIES

We are committed to **innovation that embraces intricacies and complexities.** Our exceptional tools and education programs empower physicians with the confidence they need to perform evermore challenging procedures and spark progress.

MedSphere is a system manufactured by MedSphere® International and has CE-mark. MedSphere is not available for sale in all countries. This information is provided only in respect to markets where the product is approved or cleared. MedSphere is not FDA cleared in the US for sale. MedSphere is not approved in Canada. MedSphere® International is a trademark registered in China. Please contact your Terumo local sales representative for more information.

**MEDSPHERE**  
Radiofrequency Ablation System

## EXPAND YOUR RFA OPTIONS WITH TERUMO



from ACCESS  
to CLOSURE  
INTERVENTIONAL  
ONCOLOGY



**TERUMO**  
INTERVENTIONAL  
SYSTEMS

**TERUMO**  
INTERVENTIONAL  
SYSTEMS

**Terumo Corporation**  
+81333748111

**Terumo Europe NV**  
+321 6381 211

**Terumo  
Interventional Systems EMEA**  
+331 471 609 30

EMEA Sales Offices

**Terumo Europe NV  
Africa Business Division**  
+32 16 38 13 08

**Terumo Europe NV  
Benelux Sales Division**  
Belgium:  
0800 14468  
The Netherlands:  
0800 0231938

**Terumo Europe NV  
Emerging Market Division**  
+3216 38 12 11

**Terumo Deutschland GmbH**  
+49 6196 80 230

**Terumo Deutschland GmbH  
Österreich**  
+43 2236 378020

**Terumo Deutschland GmbH  
Switzerland**  
+41 56 419 10 10

**Terumo Europe España SL**  
+34 902 10 12 98

**Terumo France S.A.S.**  
+33 130 96 13 00

**Terumo Italia S.r.l.**  
+39 06 94 80 28 00

**Terumo Russia LLC**  
+7 495 988 4740

**Terumo Sweden AB**  
+46 3174 85 880



© Registered Trademark  
Published by Terumo Europe NV

15791680.1911(0219)E CE 0197



www.terumo-europe.com





# UMBRELLA ELECTRODE

## TAKE CONTROL OF LARGE TUMOR ABLATION

- RF ablation of solid tumors using percutaneous, laparoscopic or open surgical technique
- Primary and secondary liver lesions
- Lung lesions and other solid tumors

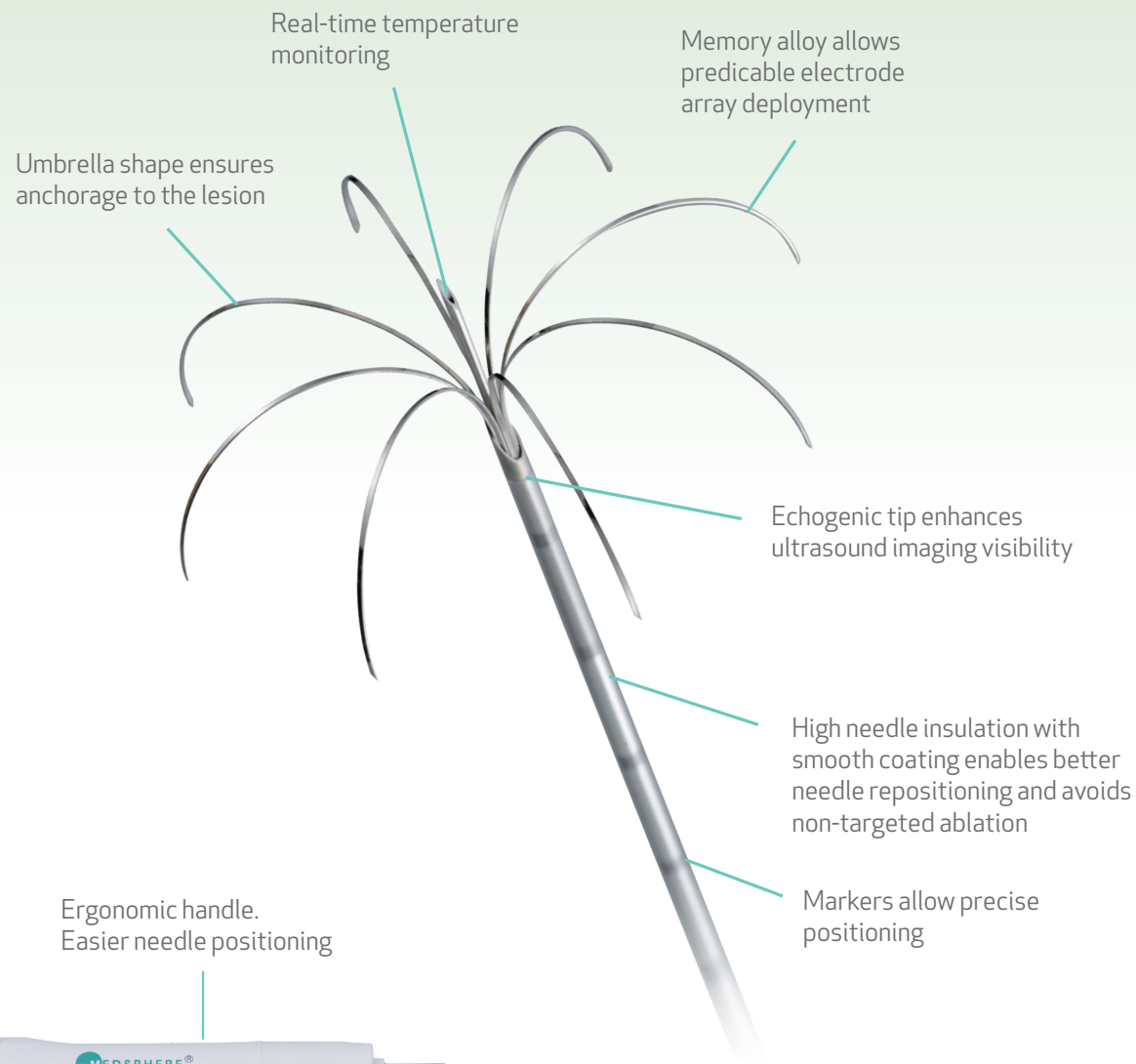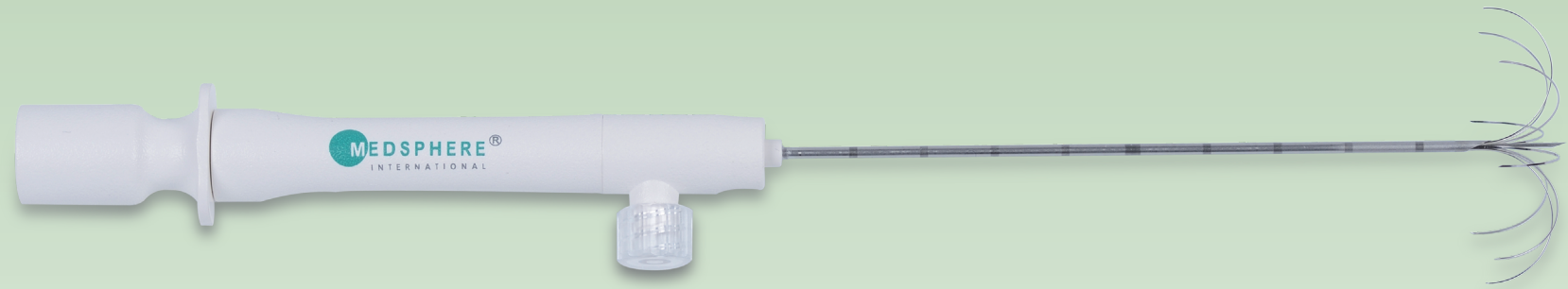


## ABLATION TREATMENT

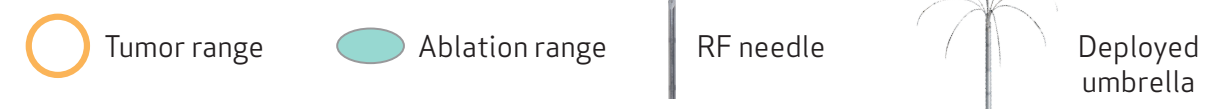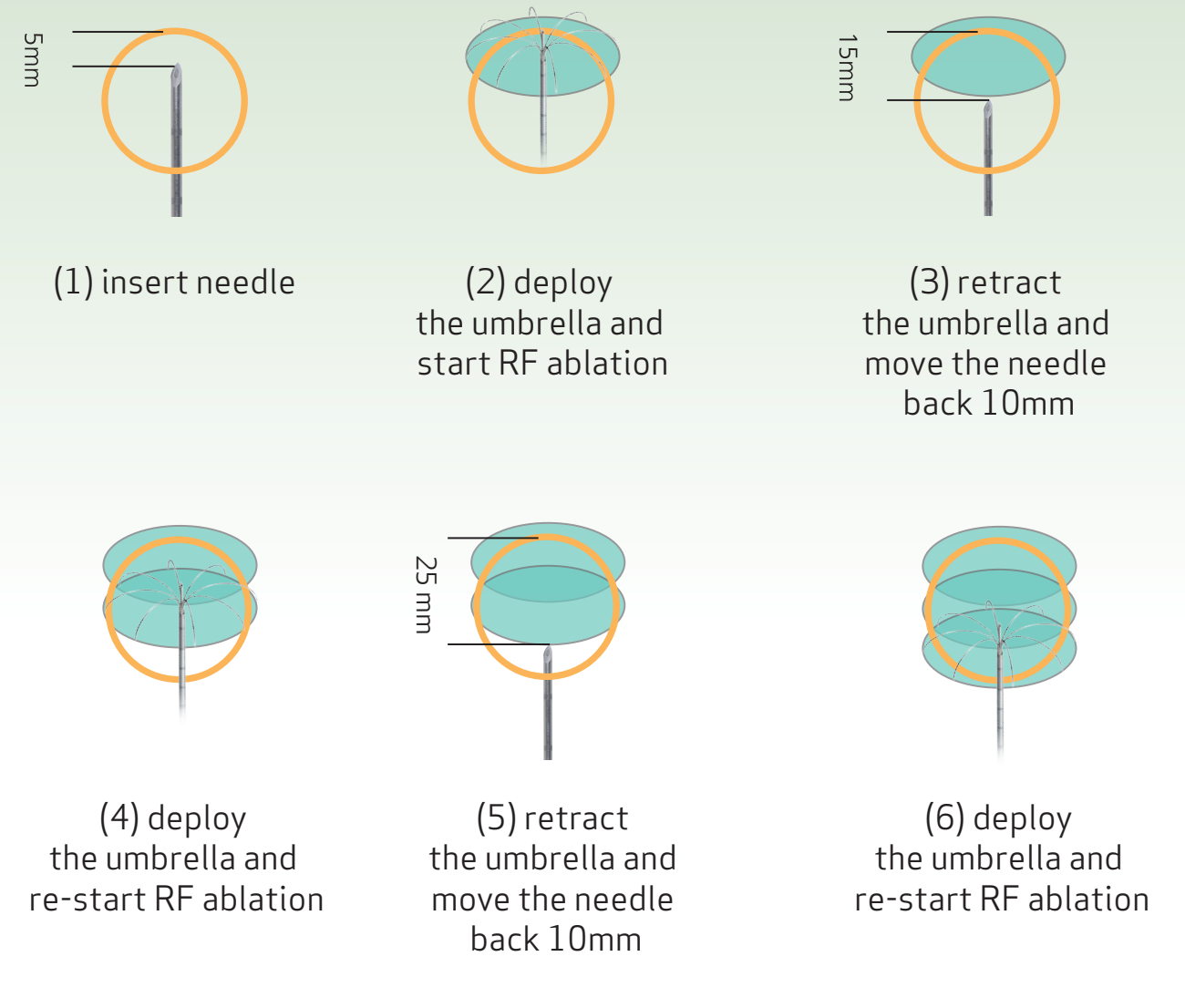

# INTERNALLY COOLED ELECTRODE

## LARGE ABLATION ZONE

- **Internally cooled tip**  
Reduced tissue adhesion and carbonization
- **19 G needle**  
The thinnest RFA needle in the market
- **Ergonomic handle**  
Easier positioning
- **High needle insulation with smooth coating**  
Enables better needle repositioning and avoids non-targeted ablation

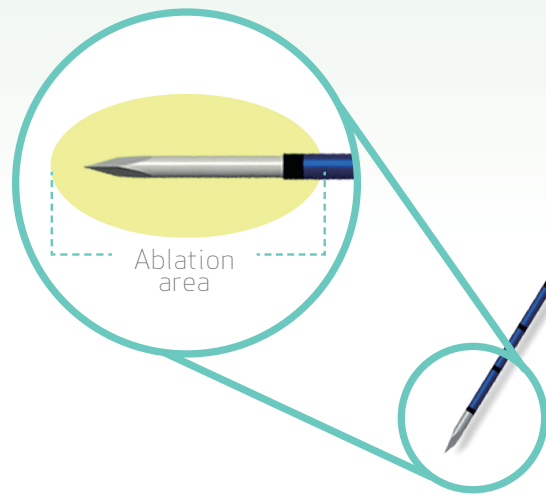

Animal Test (pig, liver) 19G, 1 cm length of tip, 30 w, 24 s. - ablation area: 14 mm\* 9 mm.




Thyroid benign nodule ablation.

# RFG-150 RF GENERATOR



## HIGH POWER GENERATOR

- **Dual mode - temperature and power control**
- **Wide power range for different lesion sizes**
- **Real-time temperature, power and impedance display**
- **Intelligent protection mechanism for enhanced safety**

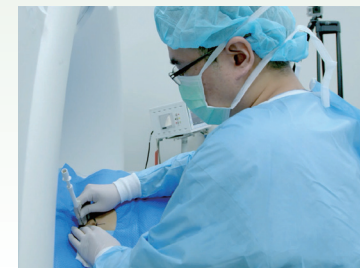

Picture from RENJI HOSPITAL



Pump



Foot Pedal

## TECHNICAL DATA


RF Power Output	5 to 150 W (1 W increment adjustments)
Operating Frequency	480 kHz
Target Temperature	50 to 120°C (1°C increment adjustments)
Measured Temperature	15 to 124°C (1°C increment adjustments)
Impedance Range	10 to 500 Ω (Measurement range 0 to 990 Ω)
Power Source	single-phase AC115-230V, 50/60 Hz, 280 VA power input
Dimensions	14cm x 32cm x 43cm (H x W x D)
Weight	7.8 KG

# ORDERING INFORMATION

## HARDWARE AND ACCESSORIES

RFG-150	RF Generator
90-090101	Peristaltic pump
1030A	Grounding pad (dispersive electrode)
90-021404	Sterile electrode connecting cable
90-080102	Tubing Set
90-020701	Non Sterile electrode connecting cable
90-030107	Dispersive electrode connecting cable
90-010205	Power cord


## UMBRELLA ELECTRODE

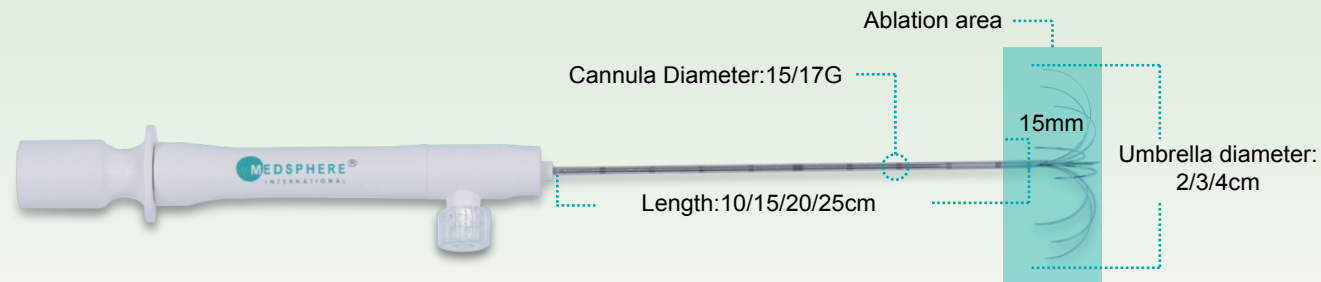

Code	Cannula Diameter	Umbrella diameter	Length	Number of tines	Thermocouple
21-843251_UPACK	15G	20mm	10cm	8	With
21-863251_UPACK	15G	20mm	15cm	8	With
21-883251_UPACK	15G	20mm	20cm	8	With
21-8A3251_UPACK	15G	20mm	25cm	8	With
21-843451_UPACK	15G	30mm	10cm	8	With
21-863451_UPACK	15G	30mm	15cm	8	With
21-883451_UPACK	15G	30mm	20cm	8	With
21-8A3451_UPACK	15G	30mm	25cm	8	With
21-843651_UPACK	15G	40mm	10cm	8	With
21-863651_UPACK	15G	40mm	15cm	8	With
21-883651_UPACK	15G	40mm	20cm	8	With
21-8A3651_UPACK	15G	40mm	25cm	8	With
21-843271_UPACK	17G	20mm	10cm	8	With
21-863271_UPACK	17G	20mm	15cm	8	With
21-883271_UPACK	17G	20mm	20cm	8	With
21-843471_UPACK	17G	30mm	10cm	8	With
21-863471_UPACK	17G	30mm	15cm	8	With
21-883471_UPACK	17G	30mm	20cm	8	With
21-843671_UPACK	17G	40mm	10cm	8	With
21-863671_UPACK	17G	40mm	15cm	8	With
21-883671_UPACK	17G	40mm	20cm	8	With

\* All kits include 2 grounding pads and 1 sterile cable.

## INTERNALLY COOLED ELECTRODE

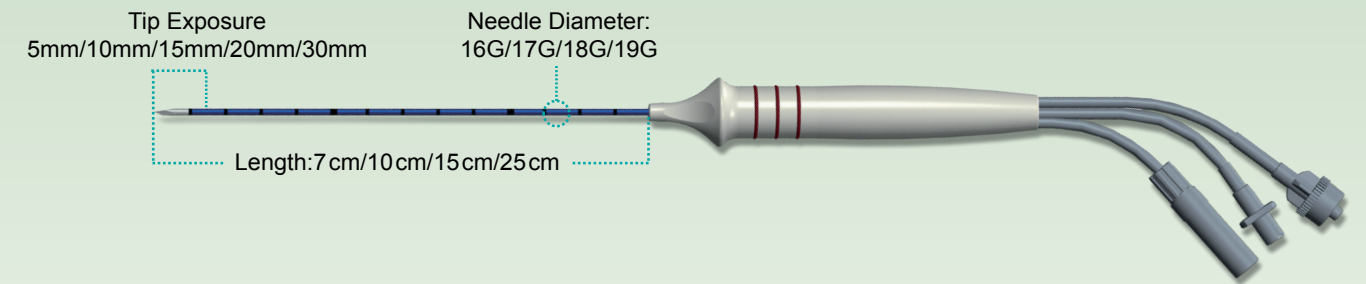

Code	Needle diameter	Tip Exposure	Length
L-342_CPACK	16G	10mm	15cm
L-442_CPACK	16G	10mm	25cm
L-343_CPACK	16G	15mm	15cm
L-344_CPACK	16G	20mm	15cm
L-444_CPACK	16G	20mm	25cm
L-445_CPACK	16G	30mm	25cm
L-332_CPACK	17G	10mm	15cm
L-333_CPACK	17G	15mm	15cm
L-334_CPACK	17G	20mm	15cm
10-181571_CPACK	17G	30mm	20cm
L-121_CPACK	18G	5mm	7cm
L-221_CPACK	18G	5mm	10cm
L-122_CPACK	18G	10mm	7cm
L-222_CPACK	18G	10mm	10cm
L-111_CPACK	19G	5mm	7cm

\*\*All kits include 2 grounding pads, 1 sterile cable and 1 tubing set.