


Renma[®]

0.035" RX PTA Balloon Dilatation Catheter



The first **R**apid **eX**change workhorse PTA Balloon

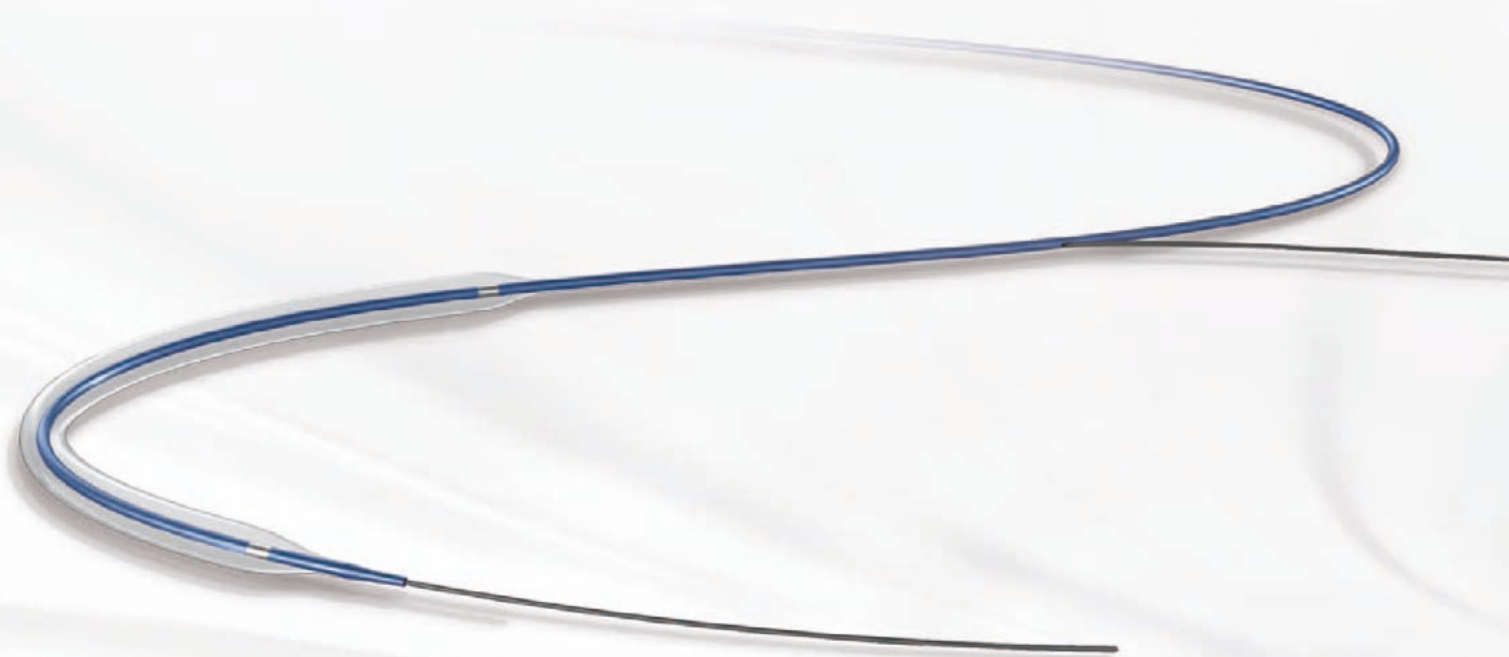

Unique 0.035" RX design

- ▶ Rapid-exchange offers a faster and more convenient exchange of devices, saves radiation dose and enables the use of short guide wires
- ▶ Perfect match with self-expanding peripheral stent Misago[®] as well in RX design

Semi-compliant 6-folded balloon with SiLX² coating

- ▶ Featuring a low crossing profile and optimal re-folding for large diameters (Ø 7-10mm)
- ▶ SiLX² coating for smooth trackability and safe retrieval after multiple inflations

Dedicated shaft design


- ▶ Combines the pushability of a stainless steel hypotube shaft with the kink resistance of a polymer shaft, even for cross-over procedures in sharp angulated bifurcations, maintains very fast deflation times

 **TERUMO[®]**

Specifications

Catheter design: Rapid Exchange, coaxial polymer/hypotube shaft
 Guide wire compatibility: 0.035" (0.89 mm)
 Usable shaft length: 135 cm
 Balloon material: Polyamide (Nylon)
 Balloon folding: 3 folded (Ø 4 – 6 mm), 6 folded (Ø 7 – 10 mm)

Coating: SiLX² coating
 Shaft diameter: Proximal 2.6Fr
 Distal 5.0-5.2Fr
 Introducer compatibility: 6Fr (2.01mm)
 Nominal Pressure (NBP): 6 atm/ 608 kPa

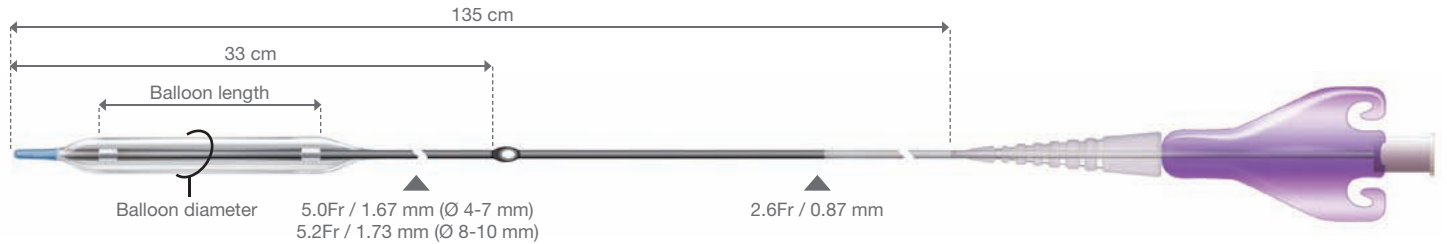

Ordering information

Renma® Order Codes

Balloon diameter	Balloon length				
	20 mm	40 mm	60 mm	80 mm	100 mm
4.0 mm	BD-C40201DN	BD-C40401DN	BD-C40601DN	BD-C40801DN	BD-C40A01DN
5.0 mm	BD-C50201DN	BD-C50401DN	BD-C50601DN	BD-C50801DN	BD-C50A01DN
6.0 mm	BD-C60201DN	BD-C60401DN	BD-C60601DN	BD-C60801DN	BD-C60A01DN
7.0 mm	BD-C70201DN	BD-C70401DN	BD-C70601DN	BD-C70801DN	BD-C70A01DN
8.0 mm	BD-C80201DN	BD-C80401DN	BD-C80601DN	-	-
9.0 mm	BD-C90201DN	BD-C90401DN	BD-C90601DN	-	-
10.0 mm	BD-CA0201DN	BD-CA0401DN	BD-CA0601DN	-	-

Rated Burst Pressure RBP

Balloon length				
20 mm	40 mm	60 mm	80 mm	100 mm
16 atm/1621 kPa			15 atm/ 1520 kPa	
14 atm/1419 kPa			13 atm/ 1317 kPa	
13 atm/1317 kPa			12 atm/ 1216 kPa	
12 atm/1216 kPa			12 atm/ 1216 kPa	
12 atm/1216 kPa			-	-
11 atm/1115 kPa			-	-
11 atm/1115 kPa			-	-



Terumo Corporation
 44-1, 2-chome, Hatagaya
 Shibuya-ku, Tokyo 151-0072,
 Japan
 Tel.: +81 3 3374 8111
 Fax: +81 3 3374 8196

European Headquarters

Terumo Europe N.V.
 Interleuvenlaan 40
 3001 Leuven, Belgium
 Tel.: +32 16 38 12 11
 Fax: +32 16 40 02 49

**Terumo Europe
 Interventional Systems
 Liaison Office**
 Immeuble Krebs
 8, rue Jean-Pierre Timbaud
 CS 10720
 78067 Saint Quentin
 en Yvelines Cedex,
 France
 Tel.: +33 139 30 95 75
 Fax: +33 139 30 95 70

European Offices

**Terumo Europe N.V.
 Benelux Sales Division**
 Interleuvenlaan 40
 3001 Leuven,
 Belgium
 Tel.: +32 16 39 25 80
 Fax: +32 16 39 25 99
 The Netherlands:
 Tel.: 0800 0220396
 Fax: 0800 0220414

Laboratoires Terumo France S.A.
 Parc Ariane - Bât. Uranus
 1, rue Hélène Boucher
 78284 Guyancourt Cedex,
 France
 Tel.: +33 130 96 13 00
 Fax: +33 130 43 60 85

Terumo Italia S.r.l.
 Via Simone Martini, 143 / 145
 00142 Roma RM, Italy
 Tel.: +39 0651 96 14 20
 Fax: +39 0650 30 407

Terumo UK Ltd
 Tamesis, The Causeway
 Egham, Surrey
 TW20 9AW, United Kingdom
 Tel.: +44 1784 476 200
 Fax: +44 1784 476 201

Terumo Deutschland GmbH
 Ludwig-Erhard-Strasse 6
 65760 Eschborn, Germany
 Tel.: +49 6196 80 230
 Fax: +49 6196 80 23 200

Terumo Deutschland GmbH
 Zweigniederlassung Spreitenbach
 Bodenackerstrasse 3
 8957 Spreitenbach, Switzerland
 Tel.: +41 56 419 10 10
 Fax: +41 56 419 10 11

Terumo Europe España SL
 Avda. Juan Carlos I, Nº 13 - 7ª Planta
 Edificio Torre La Garena
 28806 Alcalá de Henares (Madrid), Spain
 Tel.: +34 9021 01 298
 Fax: +34 9021 01 358

Terumo Sweden AB
 Nya Varvet, Byggnad 90
 Sven Källfelts Gata 18
 426 71 Västra Frölunda, Sweden
 Tel.: +46 3174 85 880
 Fax: +46 3174 85 890

**Terumo Denmark
 Filial of Terumo Sweden AB**
 Tel.: +45 7020 93 80
 Fax: +45 7020 94 80

**Terumo Europe N.V.
 Moscow Representative Office**
 BC Northern Tower
 13Fr., 10 Testovskaya Street
 Moscow 123317, Russia
 Tel.: +7 495 988 4740
 Fax: +7 495 988 4739

**Addresses subject to change,
 please consult our Website.**

© Registered Trademark.
 Published by Terumo Europe N.V.