

**FEMORAL INTRODUCER KITS**

YOUR STANDARD FOR  
**VASCULAR ACCESS**

from ACCESS  
to **CLOSURE**



 **TERUMO**  
INTERVENTIONAL  
SYSTEMS

# **RADIFOCUS™ *Introducer II***

Standard Introducer Kit





## **Cross-cut hemostasis valve**

Effectively protects against blood reflux and air aspiration

## **Total Integrated Fit (TIF) tip tapering**

Optimal tapering design at the tip of the sheath and dilator for smooth penetration

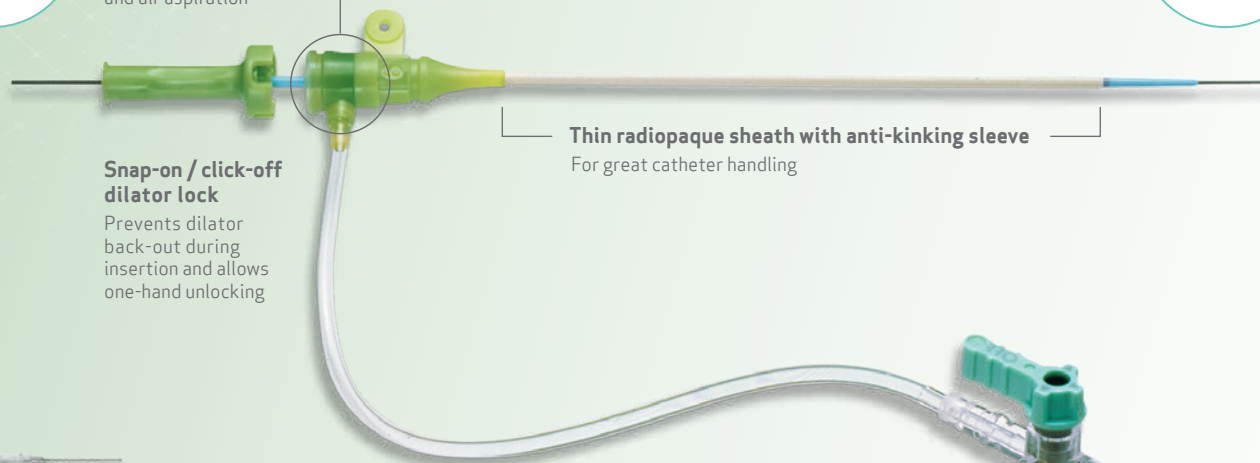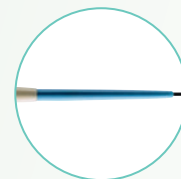

## **Snap-on / click-off dilator lock**

Prevents dilator back-out during insertion and allows one-hand unlocking

## **Thin radiopaque sheath with anti-kinking sleeve**

For great catheter handling



## **Wide variety of kit variations providing all elements for quick vessel access**

- 4 - 11 Fr sheaths
- 5 - 25 cm lengths
- Surflash or micro-puncture metal needle



## PRODUCT SPECIFICATIONS Radifocus™ Introducer II - Standard

| Kit type                                                           | Code number | Inner sheath diameter | Sheath length | Compatible guidewire (max) | Mini guidewire type and O.D.  | Entry venous cathether | Syringe | Color code                                                                            |
|--------------------------------------------------------------------|-------------|-----------------------|---------------|----------------------------|-------------------------------|------------------------|---------|---------------------------------------------------------------------------------------|
| A kit in tray with plastic guidewire + Surflash plastic IV cannula | RS*A40G10SQ | 4 Fr                  | 10 cm         | 0.025"                     | 45 cm Straight, 0.64 mm       | 20 G x 2"              | Yes     |    |
|                                                                    | RS*A40K10AQ | 4 Fr                  | 10 cm         | 0.035"                     | 45 cm Angled, 0.87 mm         | 18 G x 2 1/2"          | Yes     |    |
|                                                                    | RS*A50K10SQ | 5 Fr                  | 10 cm         | 0.035"                     | 45 cm Straight, 0.87 mm       | 18 G x 2 1/2"          | Yes     |    |
|                                                                    | RS*A50K10AQ | 5 Fr                  | 10 cm         | 0.035"                     | 45 cm Angled, 0.87 mm         | 18 G x 2 1/2"          | Yes     |    |
|                                                                    | RS*A60K10SQ | 6 Fr                  | 10 cm         | 0.035"                     | 45 cm Straight, 0.87 mm       | 18 G x 2 1/2"          | Yes     |    |
|                                                                    | RS*A60K10AQ | 6 Fr                  | 10 cm         | 0.035"                     | 45 cm Angled, 0.87 mm         | 18 G x 2 1/2"          | Yes     |    |
|                                                                    | RS*A70K10SQ | 7 Fr                  | 10 cm         | 0.035"                     | 45 cm Straight, 0.87 mm       | 18 G x 2 1/2"          | Yes     |    |
|                                                                    | RS*A80K10SQ | 8 Fr                  | 10 cm         | 0.035"                     | 45 cm Straight, 0.87 mm       | 18 G x 2 1/2"          | Yes     |    |
|                                                                    | RS*A90K10SQ | 9 Fr                  | 10 cm         | 0.035"                     | 45 cm Straight, 0.87 mm       | 18 G x 2 1/2"          | Yes     |    |
|                                                                    | RS*A10K10SQ | 10 Fr                 | 10 cm         | 0.035"                     | 45 cm Straight, 0.87 mm       | 18 G x 2 1/2"          | Yes     |    |
|                                                                    | RS*A11K10SQ | 11 Fr                 | 10 cm         | 0.035"                     | 45 cm Straight, 0.87 mm       | 18 G x 2 1/2"          | Yes     |    |
| B kit in tray with plastic guidewire                               | RS*B40G10SQ | 4 Fr                  | 10 cm         | 0.025"                     | 45 cm Straight, 0.64 mm       | No                     | No      |    |
|                                                                    | RS*B40K10AQ | 4 Fr                  | 10 cm         | 0.035"                     | 45 cm Angled, 0.87 mm         | No                     | No      |    |
|                                                                    | RS*B40K10SQ | 4 Fr                  | 10 cm         | 0.035"                     | 45 cm Straight, 0.87 mm       | No                     | No      |    |
|                                                                    | RS*B50N10AQ | 5 Fr                  | 10 cm         | 0.038"                     | 45 cm Angled, 0.97 mm         | No                     | No      |    |
|                                                                    | RS*B50N10SQ | 5 Fr                  | 10 cm         | 0.038"                     | 45 cm Straight, 0.97 mm       | No                     | No      |    |
|                                                                    | RS*B50N25AQ | 5 Fr                  | 25 cm         | 0.038"                     | 80 cm Angled, 0.97 mm         | No                     | No      |    |
|                                                                    | RS*B60N10SQ | 6 Fr                  | 10 cm         | 0.038"                     | 45 cm Straight, 0.97 mm       | No                     | No      |    |
|                                                                    | RS*B60N25AQ | 6 Fr                  | 25 cm         | 0.038"                     | 80 cm Angled, 0.97 mm         | No                     | No      |    |
|                                                                    | RS*B70N10SQ | 7 Fr                  | 10 cm         | 0.038"                     | 45 cm Straight, 0.97 mm       | No                     | No      |    |
|                                                                    | RS*B70N25AQ | 7 Fr                  | 25 cm         | 0.038"                     | 80 cm Angled, 0.97 mm         | No                     | No      |    |
|                                                                    | RS*B80N10SQ | 8 Fr                  | 10 cm         | 0.038"                     | 45 cm Straight, 0.97 mm       | No                     | No      |    |
|                                                                    | RS*B80N25AQ | 8 Fr                  | 25 cm         | 0.038"                     | 80 cm Angled, 0.97 mm         | No                     | No      |    |
|                                                                    | RS*B90N10SQ | 9 Fr                  | 10 cm         | 0.038"                     | 45 cm Straight, 0.97 mm       | No                     | No      |    |
|                                                                    | RS*B90N25AQ | 9 Fr                  | 25 cm         | 0.038"                     | 80 cm Angled, 0.97 mm         | No                     | No      |    |
|                                                                    | RS*B11N10SQ | 11 Fr                 | 10 cm         | 0.038"                     | 45 cm Straight, 0.97 mm       | No                     | No      |    |
| B kit in tray with spring double end guidewire                     | RS*B40K10MQ | 4 Fr                  | 10 cm         | 0.035"                     | 45 cm Spring, J type, 0.89 mm | No                     | No      |    |
|                                                                    | RS*B50K10MQ | 5 Fr                  | 10 cm         | 0.035"                     | 45 cm Spring, J type, 0.89 mm | No                     | No      |    |
|                                                                    | RS*B50N10MQ | 5 Fr                  | 10 cm         | 0.038"                     | 45 cm Spring, J type, 0.97 mm | No                     | No      |    |
|                                                                    | RS*B60K10MQ | 6 Fr                  | 10 cm         | 0.035"                     | 45 cm Spring, J type, 0.89 mm | No                     | No      |    |
|                                                                    | RS*B60N10MQ | 6 Fr                  | 10 cm         | 0.038"                     | 45 cm Spring, J type, 0.97 mm | No                     | No      |    |
|                                                                    | RS*B70K10MQ | 7 Fr                  | 10 cm         | 0.035"                     | 45 cm Spring, J type, 0.89 mm | No                     | No      |    |
|                                                                    | RS*B70N10MQ | 7 Fr                  | 10 cm         | 0.038"                     | 45 cm Spring, J type, 0.97 mm | No                     | No      |    |
|                                                                    | RS*B80K10MQ | 8 Fr                  | 10 cm         | 0.035"                     | 45 cm Spring, J type, 0.89 mm | No                     | No      |    |
|                                                                    | RS*B80N10MQ | 8 Fr                  | 10 cm         | 0.038"                     | 45 cm Spring, J type, 0.97 mm | No                     | No      |    |
|                                                                    | RS*B90K10MQ | 9 Fr                  | 10 cm         | 0.035"                     | 45 cm Spring, J type, 0.89 mm | No                     | No      |    |
|                                                                    | RS*B90N10MQ | 9 Fr                  | 10 cm         | 0.038"                     | 45 cm Spring, J type, 0.97 mm | No                     | No      |    |
|                                                                    | RS*B10K10MQ | 10 Fr                 | 10 cm         | 0.035"                     | 45 cm Spring, J type, 0.89 mm | No                     | No      |    |
|                                                                    | RS*B10N10MQ | 10 Fr                 | 10 cm         | 0.038"                     | 45 cm Spring, J type, 0.97 mm | No                     | No      |    |
|                                                                    | RS*B11K10MQ | 11 Fr                 | 10 cm         | 0.035"                     | 45 cm Spring, J type, 0.89 mm | No                     | No      |  |
|                                                                    | RS*B11N10MQ | 11 Fr                 | 10 cm         | 0.038"                     | 45 cm Spring, J type, 0.97 mm | No                     | No      |  |

Continued on following page

Continued from previous page

## PRODUCT SPECIFICATIONS Radifocus™ Introducer II - Standard

| Kit type                                                                | Code number  | Inner sheath diameter | Sheath length | Compatible guidewire (max) | Mini guidewire type and O.D.    | Entry venous catheter     | Syringe | Color code                                                                          |
|-------------------------------------------------------------------------|--------------|-----------------------|---------------|----------------------------|---------------------------------|---------------------------|---------|-------------------------------------------------------------------------------------|
| <b>B kit in pouch with spring double end guidewire</b>                  | RS*B40K10MR  | 4 Fr                  | 10 cm         | 0.035"                     | 45 cm Spring, J Angled, 0.89 mm | No                        | No      |  |
|                                                                         | RS*B50N10MRD | 5 Fr                  | 10 cm         | 0.038"                     | 45 cm Spring, J Angled, 0.89 mm | No                        | No      |  |
|                                                                         | RS*B60N10MRD | 6 Fr                  | 10 cm         | 0.038"                     | 45 cm Spring, J Angled, 0.89 mm | No                        | No      |  |
|                                                                         | RS*B70N10MRD | 7 Fr                  | 10 cm         | 0.038"                     | 45 cm Spring, J Angled, 0.89 mm | No                        | No      |  |
|                                                                         | RS*B80N10MRD | 8 Fr                  | 10 cm         | 0.038"                     | 45 cm Spring, J Angled, 0.89 mm | No                        | No      |  |
|                                                                         | RS*B90N10MRD | 9 Fr                  | 10 cm         | 0.038"                     | 45 cm Spring, J Angled, 0.89 mm | No                        | No      |  |
|                                                                         | RS*B10N10MRD | 10 Fr                 | 10 cm         | 0.038"                     | 45 cm Spring, J Angled, 0.89 mm | No                        | No      |  |
|                                                                         | RS*B11N10MRD | 11 Fr                 | 10 cm         | 0.038"                     | 45 cm Spring, J Angled, 0.89 mm | No                        | No      |  |
| <b>C kit in pouch without guidewire</b>                                 | RS*C40G10NR  | 4 Fr                  | 10 cm         | 0.025"                     | No mini guidewire               | No                        | No      |  |
|                                                                         | RS*C40K10NR  | 4 Fr                  | 10 cm         | 0.035"                     | No mini guidewire               | No                        | No      |  |
|                                                                         | RS*C40K25NR  | 4 Fr                  | 25 cm         | 0.035"                     | No mini guidewire               | No                        | No      |  |
|                                                                         | RS*C50N10NR  | 5 Fr                  | 10 cm         | 0.038"                     | No mini guidewire               | No                        | No      |  |
|                                                                         | RS*C60N10NR  | 6 Fr                  | 10 cm         | 0.038"                     | No mini guidewire               | No                        | No      |  |
|                                                                         | RS*C60N25NR  | 6 Fr                  | 25 cm         | 0.038"                     | No mini guidewire               | No                        | No      |  |
|                                                                         | RS*C70N10NR  | 7 Fr                  | 10 cm         | 0.038"                     | No mini guidewire               | No                        | No      |  |
|                                                                         | RS*C70N25NR  | 7 Fr                  | 25 cm         | 0.038"                     | No mini guidewire               | No                        | No      |  |
|                                                                         | RS*C80N10NR  | 8 Fr                  | 10 cm         | 0.038"                     | No mini guidewire               | No                        | No      |  |
|                                                                         | RS*C80N25NR  | 8 Fr                  | 25 cm         | 0.038"                     | No mini guidewire               | No                        | No      |  |
|                                                                         | RS*C90N10NR  | 9 Fr                  | 10 cm         | 0.038"                     | No mini guidewire               | No                        | No      |  |
|                                                                         | RS*C10N10NR  | 10 Fr                 | 10 cm         | 0.038"                     | No mini guidewire               | No                        | No      |  |
|                                                                         | RS*C11N10NR  | 11 Fr                 | 10 cm         | 0.038"                     | No mini guidewire               | No                        | No      |  |
| <b>R kit in tray with spring double end guidewire + metallic needle</b> | RS*R40K10MQ  | 4 Fr                  | 10 cm         | 0.035"                     | 45 cm Spring, J Type, 0.89 mm   | 18 G x 2 <sup>3/4</sup> " | No      |  |
|                                                                         | RS*R50N10MQ  | 5 Fr                  | 10 cm         | 0.038"                     | 45 cm Spring, J Type, 0.97 mm   | 18 G x 2 <sup>3/4</sup> " | No      |  |
|                                                                         | RS*R60N10MQ  | 6 Fr                  | 10 cm         | 0.038"                     | 45 cm Spring, J Type, 0.97 mm   | 18 G x 2 <sup>3/4</sup> " | No      |  |
|                                                                         | RS*R70N10MQ  | 7 Fr                  | 10 cm         | 0.038"                     | 45 cm Spring, J Type, 0.97 mm   | 18 G x 2 <sup>3/4</sup> " | No      |  |
|                                                                         | RS*R80N10MQ  | 8 Fr                  | 10 cm         | 0.038"                     | 45 cm Spring, J Type, 0.97 mm   | 18 G x 2 <sup>3/4</sup> " | No      |  |
|                                                                         | RS*R90N10MQ  | 9 Fr                  | 10 cm         | 0.038"                     | 45 cm Spring, J Type, 0.97 mm   | 18 G x 2 <sup>3/4</sup> " | No      |  |

**A** Kit includes: introducer sheath, dilator, mini guidewire, entry venous catheter, 2.5 cc syringe

**B** Kit includes: introducer sheath, dilator, mini guidewire

**C** Kit includes: introducer sheath, dilator

**R** Kit includes: introducer sheath, dilator, mini guidewire, 18 G x 2<sup>3/4</sup>" entry metallic needle

When the mini guidewire has an angled tip or a J tip, the kit also contains a guide inserter.

Radifocus™ Introducer II - For Pediatric Use

| Kit type                                                           | Code number | Inner sheath diameter | Sheath length | Compatible guidewire (max) | Mini guidewire type and O.D. | Entry needle size | Syringe | Color code                                                                          |
|--------------------------------------------------------------------|-------------|-----------------------|---------------|----------------------------|------------------------------|-------------------|---------|-------------------------------------------------------------------------------------|
| A kit in tray with plastic guidewire + Surflash plastic IV cannula | RS*A40G05SQ | 4 Fr                  | 5 cm          | 0.025"                     | 45 cm Straight, 0.64 mm      | 20 G x 2"         | Yes     |  |
|                                                                    | RS*A40G07SQ | 4 Fr                  | 7 cm          | 0.025"                     | 45 cm Straight, 0.64 mm      | 20 G x 2"         | Yes     |  |
|                                                                    | RS*A50G05SQ | 5 Fr                  | 5 cm          | 0.025"                     | 45 cm Straight, 0.64 mm      | 20 G x 2"         | Yes     |  |
|                                                                    | RS*A50G07SQ | 5 Fr                  | 7 cm          | 0.025"                     | 45 cm Straight, 0.64 mm      | 20 G x 2"         | Yes     |  |
|                                                                    | RS*A60G05SQ | 6 Fr                  | 5 cm          | 0.025"                     | 45 cm Straight, 0.64 mm      | 20 G x 2"         | Yes     |  |
|                                                                    | RS*A60G07SQ | 6 Fr                  | 7 cm          | 0.025"                     | 45 cm Straight, 0.64 mm      | 20 G x 2"         | Yes     |  |
| B kit in tray with plastic guidewire                               | RS*B40G07SQ | 4 Fr                  | 7 cm          | 0.025"                     | 45 cm Straight, 0.64 mm      | No                | No      |  |
|                                                                    | RS*B50G07SQ | 5 Fr                  | 7 cm          | 0.025"                     | 45 cm Straight, 0.64 mm      | No                | No      |  |
|                                                                    | RS*B60G07SQ | 6 Fr                  | 7 cm          | 0.025"                     | 45 cm Straight, 0.64 mm      | No                | No      |  |


A Kit includes: introducer sheath, dilator, mini guidewire, entry venous catheter, 2.5 cc syringe

B Kit includes: introducer sheath, dilator, mini guidewire

# **RADIFOCUS™ *Introducer II***

Introducer Kit with Radiopaque Marker

**Platinum-iridium radiopaque marker band at the sheath tip**  
Great visibility for precise positioning

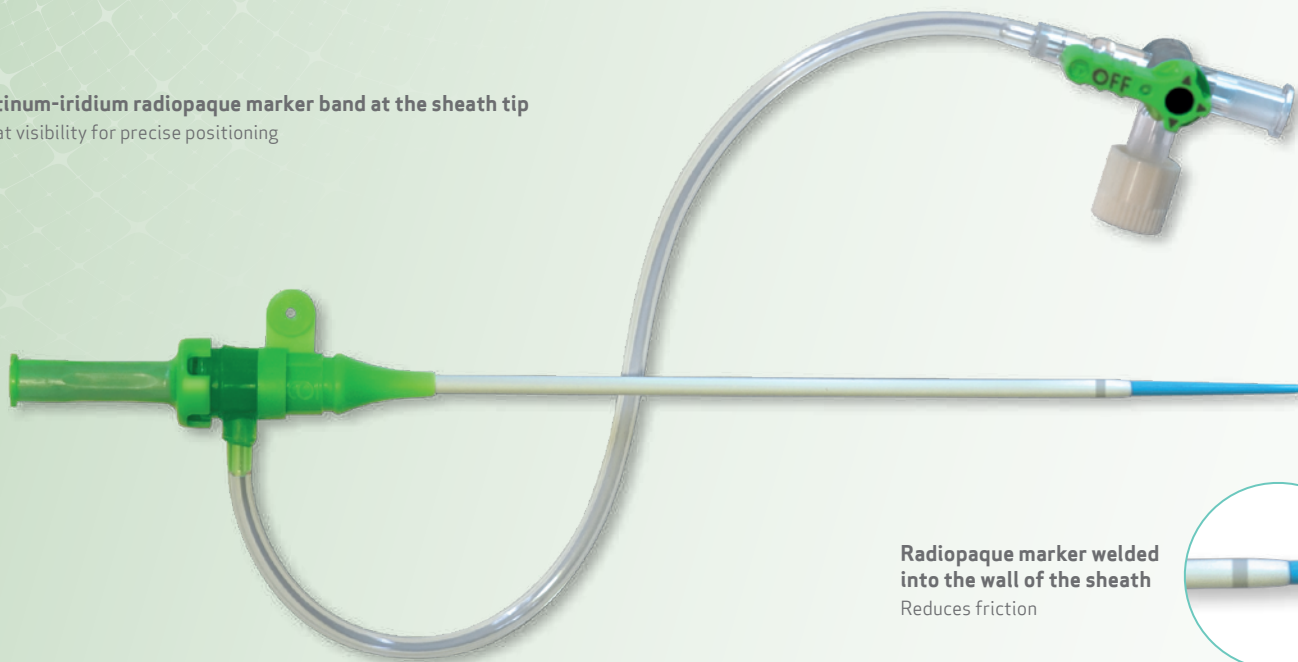

**Radiopaque marker welded  
into the wall of the sheath**  
Reduces friction

## **Radiopaque Marker Introducer Sheaths**

- available 4 - 9Fr
- lengths of 6cm, 10cm and 25cm

PRODUCT SPECIFICATIONS

B-Kit (sheath, dilator, spring mini guidewire, guide inserter)

| Code number | Inner sheath diameter | Sheath length | Compatible guidewire (max) | Mini guidewire type and O.D.    | Color code                                                                          |
|-------------|-----------------------|---------------|----------------------------|---------------------------------|-------------------------------------------------------------------------------------|
| RB*B4SK06MC | 4 Fr                  | 6 cm          | 0.035"                     | 45 cm Spring, J Angled, 0.89 mm |  |
| RB*B4SK10MC | 4 Fr                  | 10 cm         | 0.035"                     | 45 cm Spring, J Angled, 0.89 mm |  |
| RB*B5SN06MC | 5 Fr                  | 6 cm          | 0.038"                     | 45 cm Spring, J Angled, 0.97 mm |  |
| RB*B5SN10MC | 5 Fr                  | 10 cm         | 0.038"                     | 45 cm Spring, J Angled, 0.97 mm |  |
| RB*B6SN06MC | 6 Fr                  | 6 cm          | 0.038"                     | 45 cm Spring, J Angled, 0.97 mm |  |
| RB*B6SN10MC | 6 Fr                  | 10 cm         | 0.038"                     | 45 cm Spring, J Angled, 0.97 mm |  |
| RB*B7SN06MC | 7 Fr                  | 6 cm          | 0.038"                     | 45 cm Spring, J Angled, 0.97 mm |  |
| RB*B7SN10MC | 7 Fr                  | 10 cm         | 0.038"                     | 45 cm Spring, J Angled, 0.97 mm |  |
| RB*B8SN06MC | 8 Fr                  | 6 cm          | 0.038"                     | 45 cm Spring, J Angled, 0.97 mm |  |
| RB*B8SN10MC | 8 Fr                  | 10 cm         | 0.038"                     | 45 cm Spring, J Angled, 0.97 mm |  |
| RB*B9SN06MC | 9 Fr                  | 6 cm          | 0.038"                     | 45 cm Spring, J Angled, 0.97 mm |  |
| RB*B9SN10MC | 9 Fr                  | 10 cm         | 0.038"                     | 45 cm Spring, J Angled, 0.97 mm |  |

C-Kit (sheath, dilator)

| Code number | Inner sheath diameter | Sheath length | Compatible guidewire (max) | Color code                                                                            |
|-------------|-----------------------|---------------|----------------------------|---------------------------------------------------------------------------------------|
| RB*C4SK06NC | 4 Fr                  | 6 cm          | 0.035"                     |    |
| RB*C4SK10NC | 4 Fr                  | 10 cm         | 0.035"                     |    |
| RB*C4SK25NC | 4 Fr                  | 25 cm         | 0.035"                     |    |
| RB*C5SN06NC | 5 Fr                  | 6 cm          | 0.038"                     |    |
| RB*C5SN10NC | 5 Fr                  | 10 cm         | 0.038"                     |    |
| RB*C5SN25NC | 5 Fr                  | 25 cm         | 0.038"                     |    |
| RB*C6SN06NC | 6 Fr                  | 6 cm          | 0.038"                     |    |
| RB*C6SN10NC | 6 Fr                  | 10 cm         | 0.038"                     |    |
| RB*C6SN25NC | 6 Fr                  | 25 cm         | 0.038"                     |    |
| RB*C7SN06NC | 7 Fr                  | 6 cm          | 0.038"                     |    |
| RB*C7SN10NC | 7 Fr                  | 10 cm         | 0.038"                     |    |
| RB*C7SN25NC | 7 Fr                  | 25 cm         | 0.038"                     |    |
| RB*C8SN10NC | 8 Fr                  | 10 cm         | 0.038"                     |    |
| RB*C8SN25NC | 8 Fr                  | 25 cm         | 0.038"                     |   |
| RB*C9SN10NC | 9 Fr                  | 10 cm         | 0.038"                     |  |
| RB*C9SN25NC | 9 Fr                  | 25 cm         | 0.038"                     |  |



[www.terumo-europe.com](http://www.terumo-europe.com)

**Terumo Corporation**  
+81 3 3374 8111

**Terumo Europe NV**  
+32 16 38 12 11

**Terumo  
Interventional Systems EMEA**  
+33 1 47 16 09 30

*EMEA Sales Offices*

**Terumo Europe NV  
Africa Business Division**  
+32 16 38 13 08

**Terumo Europe NV  
Benelux Sales Division**  
Belgium:  
0800 14468  
The Netherlands:  
0800 0231938

**Terumo Europe NV  
Emerging Market Division**  
+32 16 38 12 11

**Terumo Deutschland GmbH**  
+49 6196 8023 0

**Terumo Deutschland GmbH  
Österreich**  
+43 2236 378020

**Terumo Deutschland GmbH  
Switzerland**  
+41 56 419 10 10

**Terumo Europe España SL**  
+34 902 10 12 98

**Terumo France S.A.S.**  
+33(0) 800 90 50 42

**Terumo Italia S.r.l.**  
+39 06 94 80 28 00

**Terumo Russia LLC**  
+7 495 988 4740

**Terumo Sweden AB**  
+46 3174 85 880

**Terumo Middle East FZE**  
+971 4 292 0200

**Terumo UK Ltd**  
+44 1276 480440

**Terumo BCT Tıbbi Cihazlar Dağıtım  
ve Hizmetleri A.Ş.**  
+90 216 645 92 00

**Terumo Poland Sp. z o.o.**  
+48 22 120 16 00

© Registered Trademark  
Published by Terumo Europe NV

IS964GB0723MVII