


# Carefully Crafted

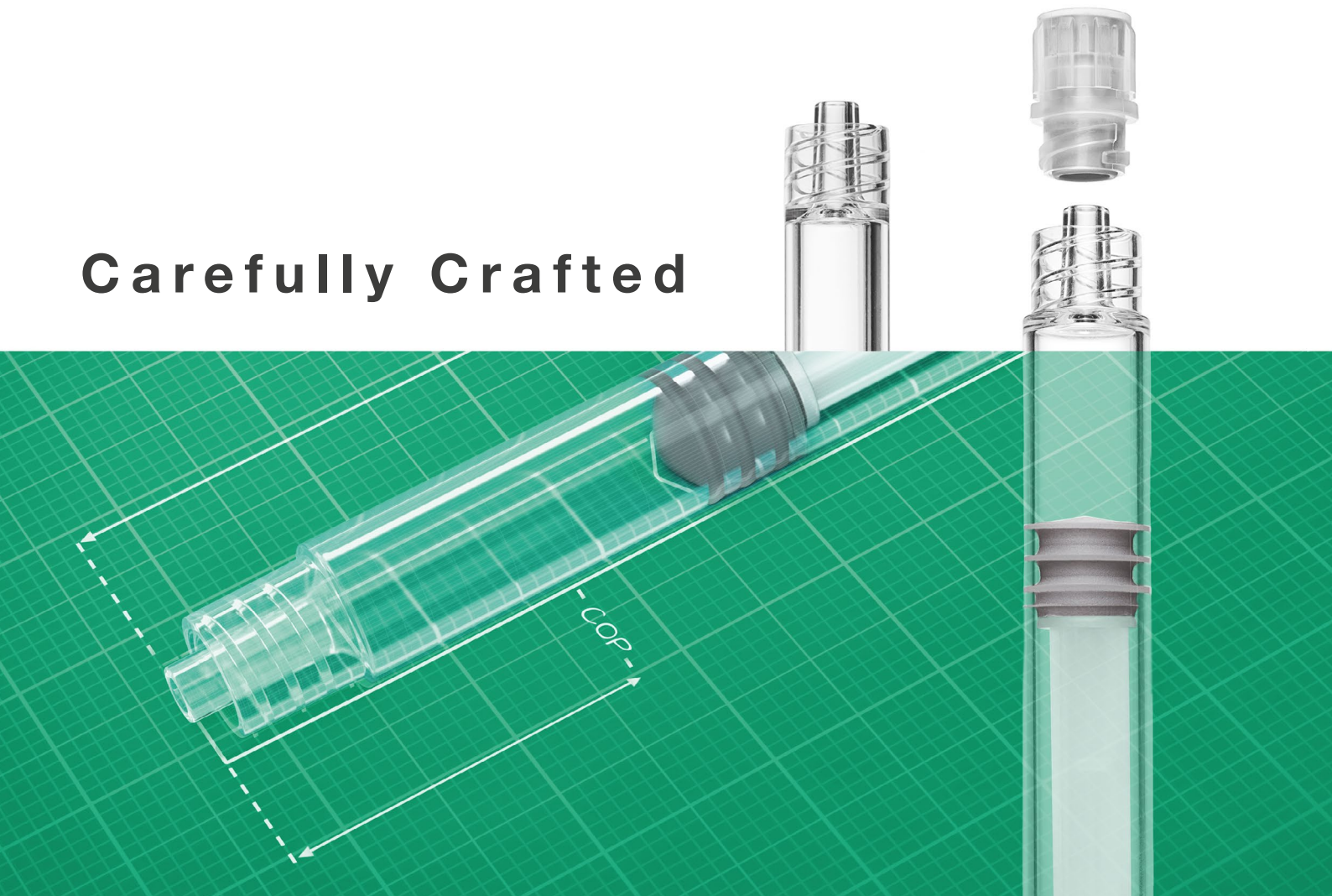


**PLAJEX**<sup>™</sup>  
Ready-to-Fill Polymer Syringe

## Polymer (COP) Pre-Fillable Syringe

### Integrated Luer Lock with Threaded Closure (TC)

- Easy grip twist off tip cap
- Luer lock for secure connection
- Luer cone/luer lock compliant with ISO 80369-7
- Large luer inner bore for connectivity with needle-free access devices
- Steam sterilization to minimize container/drug interactions

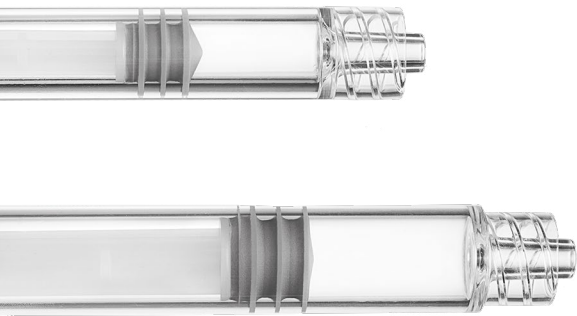

## Polymer syringe platform for reliable and effective drug delivery

### Features and benefits

- Twist-off cap for ease of use
- Both closure parts from the same elastomeric contact material (CIIR)
- Minimize extractables and leachables
- Excellent strength and clarity\*
- Dimensional accuracy
- Ready-to-fill format
- Steam sterilized



PLAJEX™ with TC Luer Lock System		
Syringe size	1 ml Long – 2.25 ml reference to ISO 11040-6	
Syringe barrel lubrication	Silicone Oil Option: 1 ml Long silicone oil-free	
Luer Lock	Integrated molded luer lock Compliant with ISO 80369-7	
Twist-off tip cap	Rubber component	Chlorobutyl rubber (CIIR)
	Rigid cap	Polypropylene
Plunger stopper	Chlorobutyl rubber (CIIR) Reference to ISO 11040-5 Optional: 1 ml Long (silicone oil-free) i-coating™	

Please note that some items may not be available in all regions, please check with your local TERUMO Pharmaceuticals Solutions contact for details.

\*data on file